



Blood Glucose Monitoring System Instructions For Use (IFU)

CONTENTS

1 Important Information About Your System	Front page	6 Testing Blood	Front page
Intended Use		Obtaining a Blood Sample	
Important Health and Safety Information		How to Test Blood	
References		System Out of Range Warning Messages	
2 System Specifications	Front page	System and Laboratory Testing	
Operating Range / Chemical Composition		7 Meter Set Up	Back page
Expected Results		Time/Date Set Up	
Symbols Chart		8 Meter Memory	Back page
3 Know Your System	Front page	View Averages (7-, 14-, 30- day)	
Meter		View Memory	
Test Strip		9 System Care	Back page
Control Solution		TRUE METRIX Control Solution	
To Attach/Remove Meter to Test Strip Vial		TRUE METRIX Blood Glucose Test Strip Care	
4 Getting Started	Front page	Meter Cleaning and Disinfecting	
Quality Control Testing		Changing Battery	
Automatic Self-Test		10 Performance Characteristics	Back page
Control Test		11 Troubleshooting	Back page
5 System Safety Information	Front page	12 Messages	Back page
Electromagnetic Compatibility			

1 IMPORTANT INFORMATION ABOUT YOUR SYSTEM

INTENDED USE
 The TRUE METRIX GO Blood Glucose Monitoring System is intended for the quantitative measurement of glucose (sugar) in fresh capillary whole blood samples drawn from the fingertip or forearm, or venous whole blood collected in only sodium heparin blood collection tubes.
 The TRUE METRIX GO System is intended for self-testing outside the body (IVD) by people with diabetes at home and for multiple-patient use in professional healthcare settings as an aid to monitor the effectiveness of diabetes control.
 The TRUE METRIX GO System should not be used for the diagnosis or screening of diabetes or for neonate (newborn) use. Alternate site (forearm) testing should be done only during steady-state times (when glucose is not changing rapidly).
 The TRUE METRIX Test Strips are for use with the TRUE METRIX GO Meter to quantitatively measure glucose (sugar) in fresh capillary whole blood samples drawn from the fingertip or forearm and venous whole blood.
 The TRUE METRIX Control Solution is for use with the TRUE METRIX GO Meter and TRUE METRIX Test Strips to check that the meter and the test strip are working together properly and that the test is performing correctly.
 The TRUE METRIX GO Meter measures the current, detects, analyzes and corrects for hematocrit and temperature, and calculates the glucose result.

Please read complete System IFU and all product Instructions for Use before using the System.

Color Codes:		
Pink - Caution: Provides information that is important for user protection and about risks for inaccurate results.	Yellow - Important: Provides important information on testing and other issues relating to testing.	Blue - Notes: Helpful hints

IMPORTANT HEALTH and SAFETY INFORMATION

- Use of the TRUE METRIX GO System in a manner not specified in this System Instructions for Use is not recommended and may affect the ability to determine true blood glucose levels.
- All meter brands perform differently. Test results from one meter brand to another may vary. This is why test results from your meter should only be compared to a lab instrument and not to another meter brand.
- Wash hands thoroughly with soap and warm water before and after handling the meter, lancing device, lancets, or test strips as contact with blood presents an infection risk.
- To help prevent false high results, wash hands before using the system to test blood, especially after fruit has been handled.
- ALL parts of the system could carry blood-borne pathogens after use, even after cleaning and disinfecting. Cleaning and disinfecting the meter and lancing device destroys most, but not necessarily all, blood-borne pathogens.
- For instructions on how to clean and disinfect the meter, see *Meter Cleaning and Disinfecting*.
- If the meter is being operated by a second person who gives testing assistance, the meter and the lancing device should be cleaned and disinfected before use by the second person. The second person should wear disposable gloves when performing testing. It is important to keep the meter and lancing device clean and disinfected.
- Alternate site (forearm) testing should not be used to calibrate continuous glucose monitors (CGMs) or used for insulin dose calculations.
- Alternate site (forearm) testing should be done only during steady-state times (when glucose is not changing rapidly).
- The System has not been tested with animals. Do not use to test blood glucose on pets.
- If there are symptoms of low or high blood glucose, check blood glucose immediately. If the result does not match how you feel, repeat the test. If the results still do not match the way you feel, call a Doctor or Healthcare Professional immediately.
 - Low blood glucose (hypoglycemia) symptoms may be trembling, sweating, intense hunger, nervousness, weakness, and trouble speaking.
 - High blood glucose (hyperglycemia) symptoms may be intense thirst, a need to urinate often, dry mouth, vomiting, and headache.
- Since any meter may fail, break, or be misplaced, always have a backup meter.
- Do not use for diagnosis of or screening for diabetes or for neonatal use.
- Inaccurate results may occur in severely hypotensive individuals or in dehydrated patients or patients in shock. Inaccurate results may occur for individuals experiencing a hyperglycemic-hyperosmolar state, with or without ketosis.

WARNING! Upon opening the test strip carton, examine the product for missing, damaged or broken parts. Ensure the test strip vial cap is securely closed. If the product is damaged or the vial cap is not closed, DO NOT use the test strips for testing; product may give inaccurate results. Use the contact information at the bottom of the page for replacement and assistance.

For best results using the TRUE METRIX GO System:

- Read all product instructions for use before testing.
- Perform a Control Test before performing a blood glucose test for the first time. Use the contact information at the bottom of the page for information on how to obtain different levels of control solution.
- Capillary whole blood from the fingertip or forearm may be used for testing with the TRUE METRIX GO System. Forearm testing should be used only during steady-state blood glucose conditions. Venous blood collected in only sodium heparin blood collection tubes may be used. Mix well before use.
- DO NOT use venous whole blood collected in sodium fluoride blood collection tubes. Blood samples containing sodium fluoride may cause false low glucose results or blood results may be read as control solution.
- Use only TRUE METRIX Test Strips and TRUE METRIX Control Solution with the TRUE METRIX GO Meter.
- Remove only one test strip at a time from the test strip vial. Recap vial immediately after removing the test strip.

NEVER reuse test strips.
 NEVER try to wipe test strips with water, alcohol, or any cleaner to remove blood or control solution to reuse test strips. Reuse of test strips will cause inaccurate results. NEVER add a second drop of sample (blood or control solution) to the test strip. Adding more sample to the test strip after testing begins gives an error message. Do not bend, cut, or alter test strips in any way.

REFERENCES

- Joslin Diabetes Center. *Goals for Blood Glucose Control* [Electronic Version]. Retrieved June 8, 2015 from <http://www.joslin.org/info/Goals-for-Blood-Glucose-Control.html>.
- FDA Public Health Notification: Use of Fingertick Devices as More than One Person Poses Risk for Transmitting Blood Borne Pathogens. *Initial Communication Update 11/29/2010* [Electronic Version]. Retrieved February 22, 2012 from <http://www.fda.gov/MedicalDevices/Safety/AlertsandNotices/ucm224025.htm>.
- U.S. Food and Drug Administration. *Blood Glucose Meters, Getting the Most Out of Your Meter*. [Electronic Version]. Retrieved July 6, 2009 from www.fda.gov/MedicalDevices/Safety/AlertsandNotices/TipsandArticlesonDevicesSafety/ucm109371.htm.
- Larsson-Cohn U. *Difference between capillary and venous blood glucose during oral glucose tolerance tests*. *Scand J Clin Lab Invest* 36:805-808, 1976.
- Data on file.
- European Committee for Standardization. *In vitro diagnostic test systems - Requirements for blood-glucose monitoring systems for self-testing in managing diabetes mellitus*. Reference number EN ISO 15197:2015(E). Brussels: European Committee for Standardization; 2015.

TRIVIDIA HEALTH
 EMERGEO EUROPE
 Westervoortdijk 60,
 6827 AT Arnhem,
 The Netherlands

CE0197 UK CA2571

TRIVIDIA HEALTH, INC.
 2400 N.W. 55TH COURT
 FORT LAUDERDALE, FL 33309 U.S.A.

UKRP / UK Importer:
Trividia Health UK Limited
 27 Old Gloucester Street
 London, United Kingdom
 WC1N 3AX
 UK: 0800 689 5035

2 SYSTEM SPECIFICATIONS

Result Range: 1.1-33.3 mmol/L	Sample: 0.5 microlitre (0.5 µL) fresh capillary whole blood from the fingertip or forearm and venous blood drawn in only sodium heparin blood collection tubes.	Test Time: Results in as little as 4 seconds	Result Value: Plasma equivalent values	Assay Method: Electrochemical	Power Supply: One 3V lithium battery #CR2032 (non-rechargeable)	Battery Life: Approximately 1,000 tests or 1 year	Automatic shut-off: After two minutes of non-use	Weight: 18 grams	Size: 4.1 cm x 3.5 cm x 2.2 cm	Memory Size: 500 results	Operating Range (Meter & Test Strips for Blood Testing): 10%-90% (Non-condensing)	Relative Humidity: 5°C-40°C	Temperature: 20%-70%	Hematocrit: Up to and including 3109 metres (10,200 feet)	Altitude: Use within specified environmental conditions only.
Chemical Composition		Test Strips: Glucose dehydrogenase-FAD (<i>Aspergillus species</i>), mediators, buffers and stabilizers.													
Control Solution: Contents: water, d-glucose, buffers, viscosity enhancing agent, salts, dye and preservatives.															
EXPECTED RESULTS															
Expected Results for people without diabetes:															
Plasma Blood Glucose Result*															
Before breakfast < 5.6 mmol/L															
Two hours after a meal < 7.8 mmol/L															
Importance of Blood Glucose Monitoring															
A Doctor or Healthcare Professional determines how often to test glucose and what the target ranges are for blood glucose results.															
Having most blood glucose results within the target range shows how well a treatment plan is working to control glucose levels. To slow or stop the complications from diabetes, keep glucose results within the target range.															
NEVER change a treatment plan without talking to a Doctor or Healthcare Professional.															

SYMBOLS:	Biological Risks
STERILIZATION:	Sterilized Using Irradiation
	Do not Resterilize
	Single Use
CONTROL:	Control
① ② ③	Control Level
SN	Serial Number
	Caution
	Use-by date
	Keep Dry
	Keep away from sunlight
	Consult Instructions for Use
	Temperature Limit
	Humidity Limitation
LOT	Lot Number
IVD	In vitro Diagnostic Medical Device
UKRP	UK Responsible Person
CE REP	Authorized Representative in the European Community
	Manufacturer
	Date of Manufacture
	Single Patient Multiple Use
	Waste Electrical and Electronic Equipment

3 KNOW YOUR SYSTEM

METER

Front of Meter

- Display - Shows test results, messages, user prompts.
- Test Port - Insert Test Strip here, with contact blocks facing up.
- Set Button - Turns meter on to view Average values and scroll through Memory, sets up date/time, adds ALT Symbol, turns meter off.
- Battery Tray - Holds battery (one non-rechargeable 3V lithium battery #CR2032).
- Meter Label - Contains serial number used to identify meter when contacting for assistance.
- Micro USB Port - Used with a cable to upload results to a computer.
- Vial Lip Cover - Locks meter onto a vial of test strips.

Back of Meter

Test Strip Vial Label

(Example only and does not reflect actual Control Test ranges)

- Lot Number (LOT) - Use for identification when contacting for assistance.
- Use-by dates (☒)
- Control Test Range - Range of numbers where Control Test result must fall to assure the system is working properly.

Control Solution Bottle Label (example)

- Lot Number (LOT) - Use for identification when contacting for assistance.
- Use-by dates (☒)
- Control Solution Level (1, 2 or 3)

TEST STRIP

Insert test strip into meter before touching Sample Tip to top of blood or control solution drop. Allow drop to be drawn into the test strip until dashes appear in the Display.

- Do not apply sample to top of test strip.
- Do not smear or scrape drop with test strip.

DO NOT insert Sample Tip into meter. This may damage meter. Do not apply more sample to the test strip after testing begins.

TO ATTACH/REMOVE METER TO TEST STRIP VIAL

To attach:

- Set test strip vial on flat surface with vial lip facing to the left.
- With Test Port facing front, place bottom of meter firmly on vial top. Meter must be seated flat on top of vial cap.
- Holding the vial, twist the meter 1/4 turn clockwise. The Test Port area on the meter should cover the vial lip if attached properly.

The meter may also be used for testing without attaching to the vial.

To remove:

- Holding the vial, twist the meter 1/4 turn counterclockwise.
- Lift off meter off the vial top.

TEST STRIP VIAL LABEL

- Lot Number (LOT) - Use for identification when contacting for assistance.
- Use-by dates (☒)
- Control Test Range - Range of numbers where Control Test result must fall to assure the system is working properly.

Write date first opened on vial label. Discard vial and unused test strips if either the open vial Use-by date or the date printed next to ☒ on vial label has passed, whichever comes first. See the test strip Instructions for Use for open vial Use-by date. Use of test strips past the Use-by dates ☒ may give incorrect test results. Discard out-of-date products and test with new products. Ranges printed on test strip vial label are for Control Test results only and are not suggested levels for blood glucose.

CONTROL SOLUTION CONTROL

- Lot Number (LOT) - Use for identification when contacting for assistance.
- Use-by dates (☒)
- Control Solution Level (1, 2 or 3)

Write date first opened on bottle label. Discard bottle and unused control solution if either 3 months after first opening or date printed next to ☒ on bottle label has passed, whichever comes first. Use of control solution past the Use-by dates ☒ may give incorrect test results. Discard out-of-date products and test with new products. Do not drink control solution.

TO ATTACH/REMOVE METER TO TEST STRIP VIAL

To attach:

- Set test strip vial on flat surface with vial lip facing to the left.
- With Test Port facing front, place bottom of meter firmly on vial top. Meter must be seated flat on top of vial cap.
- Holding the vial, twist the meter 1/4 turn clockwise. The Test Port area on the meter should cover the vial lip if attached properly.

The meter may also be used for testing without attaching to the vial.

To remove:

- Holding the vial, twist the meter 1/4 turn counterclockwise.
- Lift off meter off the vial top.

4 GETTING STARTED

The meter comes with pre-set time and date. Before using the meter for the first time or after a battery change, check time and date and update as needed.
 The meter turns on when:
 ~ a test strip is inserted into the Test Port, or
 ~ when Set Button is pressed and released (see *Meter Memory and Time/Date Set Up*).
 Meter turns off when:
 ~ the test strip is released from the meter,
 ~ the Set Button is pressed and held for 3 seconds, or
 ~ after 2 minutes of non-use.

Testing Checklist:

- Check meter for damage (cracked Display, missing button, etc.). If damage is seen, do not use meter. Use the contact information at the bottom of the page for assistance.
- Check test strip vial for damage (cracked or broken vial). Discard damaged vial and contents (test strips). Use a new vial of test strips for testing.
- Write date first opened on test strip vial label. Discard vial and unused test strips if either the open vial Use-by date or the date printed next to ☒ on vial label has passed, whichever comes first. See the test strip Instructions for Use for open vial Use-by date.
- For Control Test, make sure you have clean tissues available. A small piece of plastic wrap, aluminum foil or waxed paper may be used for control solution sample drop in the Control Test.
- Check control solution bottle for any leaks or broken cap. Discard damaged bottle and open a new one for testing.
- Write date first opened on control solution bottle label. Discard bottle if either 3 months after first opening or date printed next to ☒ has passed, whichever comes first.

QUALITY CONTROL TESTING

To assure accurate and reliable results, TRUE METRIX GO offers two kinds of quality control tests, an Automatic Self-Test and a Control Test. These tests let you know that your system is working properly and your testing technique is good.

AUTOMATIC SELF-TEST

The Automatic Self-Test lets you know if the meter and the Display are working properly. The Automatic Self-Test does not take the place of running a Control Test.

- Insert test strip into Test Port.
- Full Display appears. Check for missing segments.
- Drop Symbol begins to blink. Meter may be used for testing.

If an error message appears, the meter will not perform a test. See *Troubleshooting* or contact for assistance.

CONTROL TEST

We recommend performing Control Tests:

- before using the meter for the first time,
- for practice to ensure your testing technique is good,
- when opening a new vial of test strips,
- occasionally as a vial of test strips is used,
- if results seem unusually high or low,
- if the test strip vial has been left opened, exposed to extreme heat, cold, or humidity,
- whenever a check on the performance of the system is needed,
- if meter damage is suspected (meter was dropped, crushed, wet, etc.).

Performing a Control Test with more than one level of control solution is recommended to ensure that the system is working properly. Three levels of TRUE METRIX Control Solution are available. Use contact information at the bottom of the page for more information on how to obtain levels of control solution.

Use ONLY TRUE METRIX Control Solution for Control Test.

Ranges printed on test strip vial label being used are for Control Test results only and are not suggested levels for blood glucose. Do not drink control solution.

How To Test Control Solution

Use-by dates ☒

- Gather and check supplies. See *Getting Started - Testing Checklist*.
- Allow control solution, vial of test strips and meter to adjust to room temperature for 10 minutes.
- Wash hands. Dry thoroughly.

How To Test Control Solution, cont.

- Gently swirl or invert control solution bottle to mix. **DO NOT SHAKE.**
- Remove one test strip from vial. Close vial immediately. Use test strip quickly after taking it out of the vial.
- Insert test strip firmly into Test Port. Meter turns on.
- Keep test strip in meter until testing is finished. Do not add control solution to test strip before inserting into meter.
- Remove cap from control solution bottle. Gently squeeze a drop onto a clean tissue. Wipe off bottle tip and discard tissue. Gently squeeze a drop onto a small piece of unused aluminum foil, clear plastic wrap, or waxed paper for testing.

5 SYSTEM SAFETY INFORMATION

ELECTROMAGNETIC COMPATIBILITY

The TRUE METRIX GO meter was tested and found to comply with the electromagnetic emission and immunity requirements as specified in IEC 60601-1-2 Edition 4.1. The meter's electromagnetic emission is low. The TRUE METRIX GO has met the following requirements of IEC 60601-1-2, Edition 4.1:

EMC Test

Radiated Emissions	CISPR 11 Class B limits
Conducted Emissions Voltage	Not applicable
Radiated RF EM Fields	10v/m, 80 MHz – 2.7 GHz, 80% AM at 1 kHz
Proximity fields from RF wireless communications equipment	Per table 8.10
Power Frequency Magnetic Fields	30 A/m, 50 Hz and 60 Hz
Electrical Fast Transients / Bursts	Not applicable
Surges	Not applicable
Conducted Disturbances induced by RF fields	Not applicable
Voltage Dips and Voltage Interruptions	Not applicable
Electrostatic Discharge	+/-8kV contact; +/-15kV air discharges.

Interference from the meter to other electronically driven equipment is not anticipated. The electromagnetic environment should be evaluated prior to operation of the device. Do not use the meter in a very dry environment, especially one in which synthetic materials are present. Do not use the TRUE METRIX GO meter close to sources of strong electromagnetic radiation, as these may interfere with the proper operation of the meter. Do not use electrical equipment, including antennas, closer than 12 inches to any part of the TRUE METRIX GO meter, including cables specified by the manufacturer.

Blood Glucose Monitoring System Components

- TRUE METRIX GO Blood Glucose Monitoring Meter
- TRUE METRIX Blood Glucose Test Strips
- TRUE METRIX Control Solution
- Lancing Device
- Sterile Lancets

Kit may contain one or more of the components above. Other accessories may negatively affect EMC performance. No adverse events to the Patient and Operator are anticipated due to electromagnetic disturbances because all electrical components of the TRUE METRIX GO meter are fully enclosed.

6 TESTING BLOOD

OBTAINING A BLOOD SAMPLE

Refer to lancing device's Instructions for Use for more detailed instructions on using the lancing device.

Important Notes Regarding Forearm Testing!

- Forearm testing results cannot be used for continuous glucose meter calibration or for insulin dose calculations.
- Check with a Doctor or Healthcare Professional to see if forearm testing is appropriate.
- Results from forearm are not always the same as results from fingertip.
- Some lancing devices include a special end cap for alternate site (forearm) testing. Check lancing device Instructions for Use.
- Use fingertip instead of forearm for more accurate results:
 - Within 2 hours of eating, exercise, or taking insulin,
 - If blood sugar may be rising or falling rapidly,
 - If routine results are often fluctuating,
 - If the patient is ill or under stress,
 - If forearm results do not match how you feel,
 - If blood sugar may be low or high,
 - If symptoms of low or high blood sugar are not evident.

Wash hands thoroughly with soap and warm water before and after handling the meter, lancing device, lancets or test strips.

If the meter or lancing device is being operated by a second person who gives testing assistance, the meter and the lancing device should be cleaned and disinfected before use by another person.

For instructions on how to clean and disinfect the meter, see *Meter Cleaning and Disinfecting*.

ALL parts of the system could carry blood-borne pathogens after use, even after cleaning and disinfecting. Cleaning and disinfecting the meter and lancing device destroys most, but not necessarily all, blood-borne pathogens.

- Do not reuse lancets.
- Reuse of devices labeled for single-use may result in product contamination and infection.
- Do not use lancet if protective lancet cover is already removed or broken.
- Used test strips and lancets are considered biohazardous.
- Dispose used test strips and lancets carefully into an appropriate waste container.

To help prevent false high results, wash hands before using the system to test blood, especially after fruit has been handled.

From Fingertip

- Place end of lancing device equipped with lancet against fingertip. Lance fingertip.

From Forearm

- Rub area vigorously or apply a warm dry compress to increase blood flow.


- Set lancing device aside. To help blood drop form, lower hand to waist level and gently massage from palm to fingertip.
- Place end of lancing device equipped with a lancet firmly against forearm. Press trigger button. Apply firm pressure on lancing device for 10 seconds.

then

Allow blood drop to form before attempting to apply to the test strip.

Obtaining a Blood Sample, cont.

- After testing, recap and remove used lancet from lancing device. Discard used lancet into an appropriate waste container.

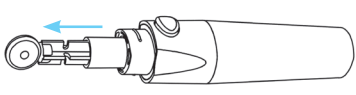


⚠ The used lancet may be biohazardous. ♻ Please discard it carefully into an appropriate waste container.

HOW TO TEST BLOOD

- Check Meter time and date before each glucose test. With meter off, press and hold the "M" Button until the full Display is shown and begins to blink. Release "M" Button. If time and date are incorrect, see *Time/Date Set Up* to set the correct time and date.

Test Strip Label

Use-by dates

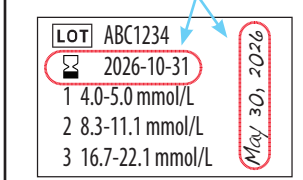


- Check supplies (see *Getting Started - Testing Checklist*).
- Allow vial of test strips and meter to adjust to room temperature for 10 minutes.




- Clean and disinfect fingertip (or forearm for alternative site testing) (see *Obtaining a Blood Sample*).
- Remove one test strip from vial. Close test strip vial immediately. Use test strip quickly after taking it out of the vial.



- Wait until Drop Symbol appears in the Display.

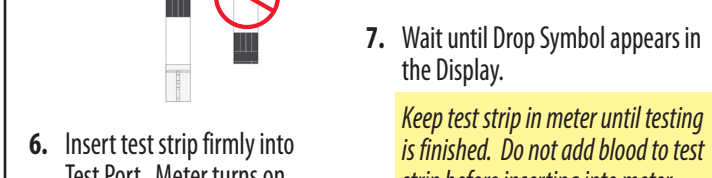


- Insert test strip firmly into Test Port. Meter turns on.

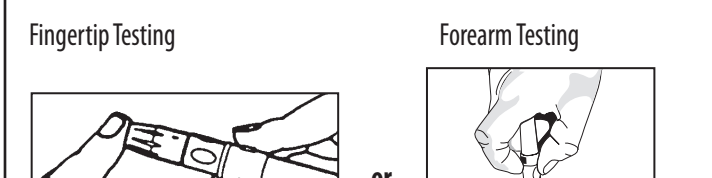

- Lance fingertip (or forearm). Allow blood drop to form (see *Obtaining a Blood Sample*).




- Remove test strip from drop when dashes appear across the meter Display. Meter is testing.

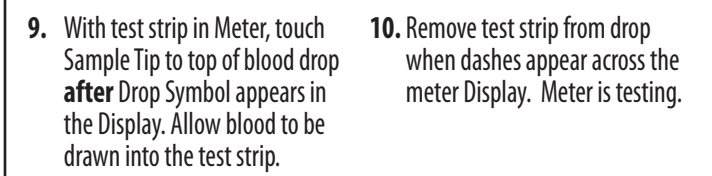

- With test strip in Meter, touch Sample Tip to top of blood drop after Drop Symbol appears in the Display. Allow blood to be drawn into the test strip.



- Remove test strip from meter and discard into an appropriate waste container. Meter turns off. Result is stored in the Memory.

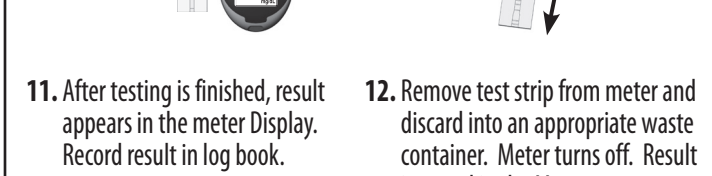

- After testing is finished, result appears in the meter Display. Record result in log book.
- If you wish to mark the result as alternate site, press the Set Button before the meter turns off. -A- appears in the meter Display.


How to Test Blood, cont.

⚠ Used test strips and lancets may be considered biohazardous. Please discard them carefully into an appropriate waste container.

- If test strip has been out of the vial too long before testing, an error message appears. Remove and discard old test strip. Use new test strip for testing.
- Do not place blood drop on top of test strip.
- Removing the test strip before result is displayed cancels the test. An error message appears. Result is not stored in Memory. Retest with a new test strip. Do not remove before result is displayed.

SYSTEM OUT OF RANGE WARNING MESSAGES

Meter reads blood glucose levels from 1.1-33.3 mmol/L. If blood test result is less than 1.1 mmol/L, "Lo" appears in meter Display. If blood test result is greater than 33.3 mmol/L, "Hi" appears in meter Display.

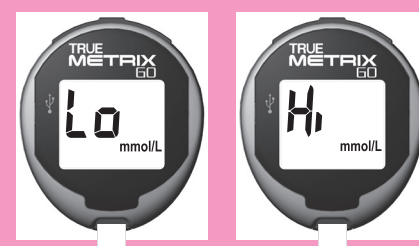

ALWAYS repeat test to confirm Low ("Lo") and High ("Hi") results. If results still display "Lo" or "Hi", call a Doctor or Healthcare Professional immediately.

~ "Lo" results are included in the Average as 1.1 mmol/L.
~ "Hi" results are included as 33.3 mmol/L.

TRUE METRIX GO SYSTEM AND LABORATORY TESTING

When comparing results between TRUE METRIX GO System and a laboratory system, TRUE METRIX GO blood tests should be performed within 30 minutes of a laboratory test. If you have recently eaten, fingertick results from the TRUE METRIX GO System can be up to 3.9 mmol/L higher than venous laboratory results.*

7 METER SET UP

TIME/DATE SET UP
If meter turns off at any time during Set Up, go back to Step #1 and begin again.



- Start with the meter off. Press and hold the Set Button until the full Display appears and begins to blink. Release Set Button.
- The time appears and the hour begins flashing. Change the number by pressing the Set Button until the desired number appears.

Pressing and releasing the Set Button only makes the numbers increase by one. Once the number reaches its limit, it resets to the lowest number. Pressing and holding the Set Button scrolls the numbers. Release Set Button when desired number is reached.



- After the correct hour appears, the number flashes for about 10 seconds before going to the minutes.
- Repeat steps 2-3 to set up minutes, month, day and year.

8 METER MEMORY

VIEW AVERAGES (7-, 14-, 30-DAY)
Averages allows you to view the average of all blood glucose results within a 7-, 14-, or 30-day period.



- With meter off, press and release Set Button.
- Display scrolls through the 7-, 14-, and 30-day Averages. Meter turns off after 2 minutes if Set Button is not pressed.

If there are no average values, three dashes are displayed for 7-, 14-, and 30-day Averages.

VIEW MEMORY

Memory stores 500 results which are displayed from most recent to oldest. When Memory is full, the oldest result is replaced with the newest result.



- After meter displays the Averages, press and release Set Button again.
- The most recent result is shown with the Memory Symbol.



- Continue to press and release the Set Button to scroll through results.



~ Only the last Control Test is saved in Memory. The Control Test result will show the Control Symbol in the lower left corner of the Display.

~ Alternate Site Blood Test results are shown with the Alternate Site Symbol in the lower left corner of the Display.

9 SYSTEM CARE

- Store system (meter, control solution, test strips) in carrying case to protect from liquids, dust and dirt.
- Store in a dry place at room temperature 4°C-30°C and at 10%-80% relative humidity (Non-condensing). **DO NOT FREEZE.**
- Allow system to sit at room temperature for 10 minutes before testing.

TRUE METRIX CONTROL SOLUTION CARE

- Write date first opened on control solution bottle label. Discard bottle and unused control solution if either 3 months after first opening or date printed next to M on label has passed, whichever comes first.
- After each use, wipe bottle tip clean and recap tightly.
- Store at room temperature 2°C-30°C. **DO NOT FREEZE.**

TRUE METRIX BLOOD GLUCOSE TEST STRIP CARE

- Store test strips in original vial only. Do not transfer test strips to new vial or store test strips outside of vial.
- Write date first opened on test strip vial label. Discard vial and unused test strips if either the open vial Use-by date or the date printed next to M on vial label has passed, whichever comes first. See the test strip Instructions for Use for open vial Use-by date. Use of test strips past the Use-by dates may give incorrect results.
- Close vial immediately after removing test strip.
- Store in a dry place at room temperature 4°C-30°C and at 10%-80% relative humidity (Non-condensing). **DO NOT FREEZE.**
- Do not store in bathroom or Do not expose to extreme heat or cold, direct sunlight or high humidity for any length of time.

METER CLEANING AND DISINFECTING

Cleaning removes blood and soil from the meter.

⚠ If the meter is being operated by a second person who provides testing assistance, the meter and lancing device should be cleaned and disinfected prior to use by the second person.

~ Do not clean the meter during a test.

How to Clean the Meter

- Wash hands thoroughly with soap and water or wear disposable gloves.




- To Clean: Make sure that the meter is off and a test strip is not inserted. Remove meter from test strip vial.

Meter may remain attached to the test strip vial during cleaning. Ensure that vial cap is completely closed before attempting to clean the system.

- Wipe meter with a clean, lint-free cloth dampened with 70% isopropyl alcohol.
- Let meter air dry thoroughly before using to test.
- Do not use bleach to clean the meter. For assistance, use the contact information at the bottom of the page.

~ Make sure no liquids enter the Test Port or any other opening in the meter.
~ If the meter remains attached to the test strip vial during cleaning, make sure that the vial cap is closed before cleaning.


How to Clean the Meter, cont.

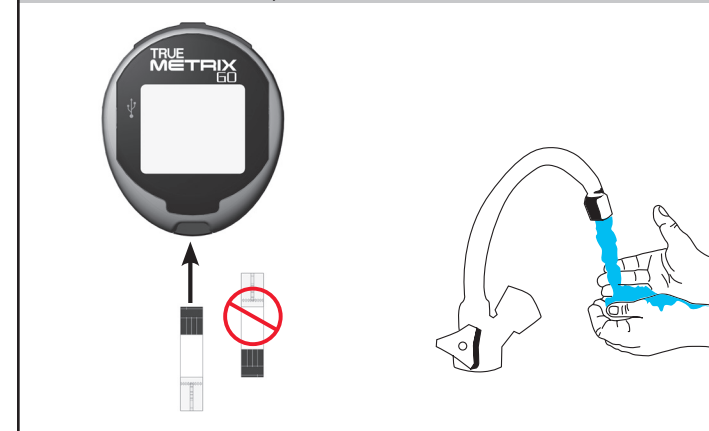

- Make sure that the system is working properly by performing an Automatic Self-Test. See *Automatic Self-Test* under *Getting Started*.
- Wash hands thoroughly after cleaning the meter.

⚠ Do not use meter and contact for assistance if:
~ Meter Display appears cloudy or any display segments are missing,
~ Markings on meter, including back meter label, are coming off or missing,
~ Set Button is hard to push on meter or does not work (see Meter Memory),
~ Unable to insert test strip into Test Port.
~ If Automatic Self-Test gives an error message.

HOW TO DISINFECT THE METER

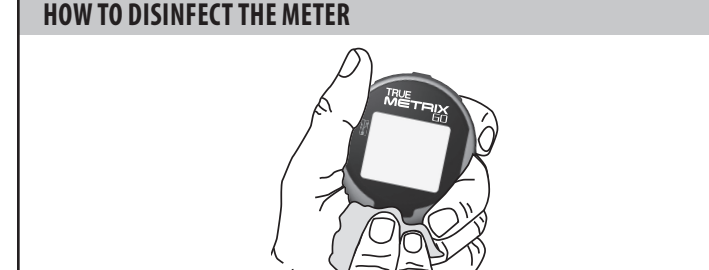

- After cleaning the meter, use a fresh wipe to wet all the outside surfaces of the meter.


~ Make sure no liquids enter the Test Port or any other opening in the meter.
~ Do not spray meter with any cleaning agents.

⚠ If the meter remains attached to the test strip vial during cleaning, make sure that the vial cap is completely closed before disinfecting.



- Keep the meter wet for 2 minutes by rewiping all of the surfaces. Fresh wipes may be used if needed to keep the meter wet. Discard used wipes.

⚠ Do not wrap or scrub the meter or the test strip vial with the wipe to keep the surfaces wet. Wrapping wipe around test strip vial or scrubbing can cause the label to become unreadable.

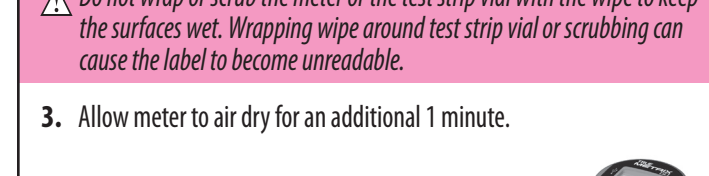

- Allow meter to air dry for an additional 1 minute.

- Reattach meter to the new test strip vial by setting the test strip vial on flat surface with vial lip facing to the left.

~ With Test Port facing front, place bottom of meter firmly on vial top. Meter must be seated flat on top of vial cap.

~ Holding the vial, twist the meter 1/4 turn clockwise. The Test Port area on the meter should cover the vial lip if attached properly.

- Wash hands thoroughly after cleaning/disinfecting the meter. Make sure your hands are completely dry.

⚠ If hands are not dry before removing a test strip from the vial, the vial may become contaminated. Test strips that have become contaminated may result in incorrect results or error messages.


LANCING DEVICE CLEANING AND DISINFECTING

Cleaning and disinfecting the lancing device is a two-step process. The first step is **Cleaning**. Cleaning removes any blood, liquids, dust or dirt. The second step is **Disinfecting**. Disinfecting removes most, but not all possible infectious agents (bacteria or viruses) including blood-borne pathogens. **Cleaning must be done before Disinfecting.**

⚠ Clean and disinfect the lancing device at least once a week. Lancing device may be cleaned and disinfected once a day for up to 3 years.

~ If blood gets on the lancing device, clean and disinfect it before anyone else touches it.

~ If someone is helping you test your blood or performs a test on you, do the cleaning and disinfecting before allowing anyone else to handle it.

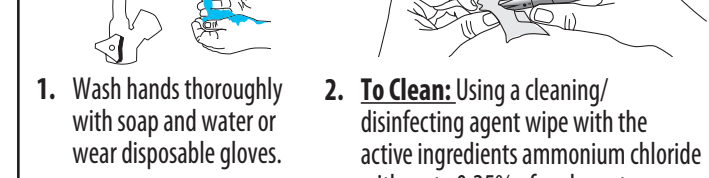

- Wash hands thoroughly with soap and water or wear disposable gloves.
- To Clean: Using a cleaning/disinfecting agent wipe with the active ingredients ammonium chloride with up to 0.25% of each quaternary ammonium compound and isopropyl alcohol (up to 55%), thoroughly wipe the entire outside surface of the lancing device with one wipe by wiping in circles 3 times. Discard the used wipe.

Lancing Device Cleaning and Disinfecting, cont.

- Repeat step 2 if the lancing device is not visually clean, or safely dispose of it according to your local regulations.
- To Disinfect: Remove a cleaning/disinfecting agent wipe with the active ingredients ammonium chloride with up to 0.25% of each quaternary ammonium compound and isopropyl alcohol (up to 55%) and thoroughly wipe the entire outside surface of the lancing device by wiping in circles 3 times. Make sure that you go over the entire lancing device.

- Using fresh wipes, if needed, make sure that all outside surfaces of the lancing device stay wet for 2 minutes by patting with the wipe. Discard used wipe. Let lancing device air dry completely before using.
- Gently pull back Arming Barrel and press the Trigger Button. You will hear a click if the lancing device is working properly.
- Wash hands thoroughly again.


⚠ Other disinfectants have not been tested. The effect of other disinfectants used interchangeably has not been tested with the lancing device. Use of disinfectants other than cleaning/disinfecting agent wipes with the active ingredients ammonium chloride with up to 0.25% of each quaternary ammonium compound and isopropyl alcohol (up to 55%) may damage lancing device.

Do not use lancing device and contact for assistance, if:
~ Markings on lancing device are coming off,
~ Trigger Button is hard to push,
~ End Cap does not go back on,
~ Arming Barrel does not click when gently pulled back.


⚠ Cleaning/disinfecting agent wipes with the active ingredients ammonium chloride with up to 0.25% of each quaternary ammonium compound and isopropyl alcohol (up to 55%) have been tested on the lancing device for a total of 1,095 cleaning and disinfecting cycles, which is equal to cleaning and disinfecting the lancing device once a day for a 3-year period. Lancing device may be cleaned and disinfected up to 7 cycles per week for 3 years.

CHANGING BATTERY

- Remove meter from top of test strip vial by holding the vial and twisting the meter 1/4 turn counterclockwise. Lift meter from vial top.
- Turn meter over until the meter label is facing up. Pull Battery Tray out until battery is exposed.

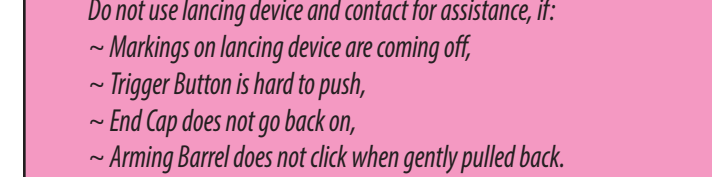

- Holding the Battery Tray over your hand, press on edge of battery until battery drops out.
- Insert new battery into Battery Tray with "+" side facing up. Slide Battery Tray back into meter.

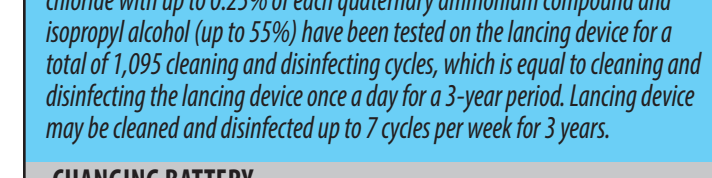

- Turn meter back over and press Set Button to turn meter on.

If meter does not turn on, open Battery Tray and check that the battery was inserted with the "+" side facing up. Close Battery Tray and repeat Step 5. If meter still does not turn on, use the contact information at the bottom of the page for assistance.

⚠ Battery may explode if mishandled. Do not dispose of battery in fire. Do not take apart or attempt to recharge battery. Dispose according to local regulations.

10 PERFORMANCE CHARACTERISTICS*

PRECISION: Precision describes the variation between results. There are two types of precision results measured - repeatability (using blood) and intermediate precision (using control solution).

Repeatability: N=100
Mean (mmol/L) 1.3 2.1 4.1 7.7 11.4 16.4 27.6
SD (mmol/L) 0.05 0.08 0.13 0.25 0.38 0.53 0.75
%CV 4.2 3.8 3.2 3.3 3.3 3.2 2.7

Intermediate Precision: N=100
Mean (mmol/L) 2.1 6.4 17.7
SD (mmol/L) 0.1 0.2 0.6
%CV 4.2 3.4 3.3

SYSTEM ACCURACY: Diabetes experts have suggested that glucose meters should agree within ±0.83 mmol/L of the medical laboratory values at glucose concentrations below 5.55 mmol/L, and within ±15% of the medical laboratory values at glucose concentrations at or above 5.55 mmol/L. The tables below show how often healthcare professionals (HCP) and users achieve these goals using capillary fingertip, capillary forearm, and venous blood samples when glucose results are not fluctuating. The laboratory reference instrument is the Yellow Springs Instrument (YSI).

Performance Characteristics*, cont.

FOR HEALTHCARE PROFESSIONALS
99.5% of TRUE METRIX GO fingertip values performed by healthcare professionals (HCP) fell within ±0.83 mmol/L of the YSI results at glucose levels <5.55 mmol/L and within ±15% at glucose levels ≥5.55 mmol/L.

Fingertip Capillary Samples (HCP vs. YSI) for glucose concentrations <5.55 mmol/L

Within ±0.28 mmol/L	Within ±0.56 mmol/L	Within ±0.83 mmol/L
94 / 156 (60.3%)	146 / 156 (93.6%)	155 / 156 (99.4%)

Fingertip Samples (HCP vs. YSI) for glucose concentrations ≥5.55 mmol/L

Within ± 5%	Within ± 10%	Within ± 15%
227 / 444 (51.1%)	383 / 444 (86.3%)	442 / 444 (99.5%)

Fingertip Samples for glucose concentrations between 1.1-33.3 mmol/L

Within ±0.83 mmol/L
597/600 (99.5%)

Parques Error Grid: 100% of individual fingertip glucose measured values performed by healthcare professionals fell within Zone A of the Parques Error Grid (PEG).

98.2% of TRUE METRIX GO forearm values performed by healthcare professionals (HCP) fell within ±0.83 mmol/L of the YSI results at glucose levels <5.55 mmol/L and within ±15% at glucose levels ≥5.55 mmol/L.

Forearm Capillary Samples (HCP vs. YSI) for glucose concentrations <5.55 mmol/L

Within ±0.28 mmol/L	Within ±0.56 mmol/L	Within ±0.83 mmol/L
28 / 62 (45.2%)	53 / 62 (85.5%)	60 / 62 (96.8%)

Forearm Samples for glucose concentrations between 1.1-33.3 mmol/L

Within ± 5%	Within ± 10%	Within ± 15%
74 / 156 (47.4%)	132 / 156 (84.6%)	154 / 156 (98.7%)

Forearm Samples for glucose concentrations between 1.1-33.3 mmol/L

Within ±0.83 mmol/L (15 mg/dL) and ±15%
214 / 218 (98.2%)

Parques Error Grid: 99.1% of individual forearm glucose measured values performed by healthcare professionals fell within Zone A and 0.9% in Zone B of the Parques Error Grid (PEG).

Venous Blood
99.1% of TRUE METRIX GO venous values performed by healthcare professionals (HCP) fell within ±0.83 mmol/L of the YSI results at glucose levels <5.55 mmol/L and within ±15% at glucose levels ≥5.55 mmol/L.

Venous Samples (HCP vs. YSI) for glucose concentrations <5.55 mmol/L

Within ±0.28 mmol/L (5 mg/dL)	Within ±0.56 mmol/L	Within ±0.83 mmol/L
61 / 90 (67.8%)	85 / 90 (94.4%)	90 / 90 (100%)

Venous Samples (HCP vs. YSI) for glucose concentrations ≥5.55 mmol/L

Within ± 5%	Within ± 10%	Within ± 15%
66 / 130 (50.8%)	122 / 130 (93.8%)	128 / 130 (98.5%)

Venous Samples for glucose concentrations between 1.1-33.3 mmol/L

Within ±0.83 mmol/L (15 mg/dL) and ±15%
218/220 (99.1%)

Parques Error Grid: 100% of individual venous glucose measured values performed by healthcare professionals fell within Zone A of the Parques Error Grid (PEG).

FOR CONSUMERS
99% of TRUE METRIX GO fingertip values performed by users fell within ±0.83 mmol/L (15 mg/dL) of the YSI results at glucose levels <5.55 mmol/L and within ±15% at glucose levels ≥5.55 mmol/L.

Fingertip Samples (User vs. YSI) for glucose concentrations <5.55 mmol/L

Within ±0.28 mmol/L	Within ±0.56 mmol/L	Within ±0.83 mmol/L
13 / 17 (76.5%)	17 / 17 (100%)	17/17 (100%)

Fingertip Samples (User vs. YSI) for glucose concentrations ≥5.55 mmol/L

Within ± 5%	Within ± 10%	Within ± 15%
46 / 83 (55.4%)	73 / 83 (88.0%)	82/83 (98.8%)

Parques Error Grid: 100% of individual fingertip glucose measured values performed by users fell within Zone A of the Parques Error Grid (PEG).

98.2% of TRUE METRIX GO forearm values performed by users fell within ±0.83 mmol/L of the YSI results at glucose levels <5.55 mmol/L and within ±15% at glucose levels ≥5.55 mmol/L.

Forearm Samples (User vs. YSI) for glucose concentrations <5.55 mmol/L

Within ±0.28 mmol/L	Within ±0.56 mmol/L	Within ±0.83 mmol/L
13 / 31 (41.9%)	22 / 31 (71.0%)	31/31 (100%)

Forearm Samples (User vs. YSI) for glucose concentrations ≥5.55 mmol/L

Within ± 5%	Within ± 10%	Within ± 15%
34 / 78 (43.6%)	64 / 78 (82.1%)	76 / 78 (97.4%)

Forearm Samples for glucose concentrations between 1.1-33.3 mmol/L

Within ±0.83 mmol/L
99/100 (99.0%)

Parques Error Grid: 100% of individual forearm glucose measured values performed by users fell within Zone A of the Parques Error Grid (PEG).